

53A SERIES

Single Stack Filter Housing

The **53A Series** single cartridge housing is ideal for use in sensitive gas, steam and pulsating pressure environments. With the ability to accept both standard Double Open End and Single Open End cartridges, this single stack housing uses a DIN threaded closure that has proved successful in a wide variety of applications such as aerosol filling lines and nitrogen purge systems for beverage and pharmaceutical applications.

The **53A Series** will accept the full range of Amazon and industry standard cartridges ranging from 0.1 to 500µm with nominal lengths of 5, 10, 20 and 30".




Double Open Ended

Single Open Ended

Product Features

- 3 piece construction with loose ring
- Suitable for gas, steam, and pulsating pressure environments
- Lightweight, high specification investment cast head and top ring
- Easy to use tie rod assembly design ensures cartridge centralisation (for a reliable knife edge seal) that enables the use of either 9.75" or 10" cartridge length configurations or multiples) in the same housing.
- Dual purpose head accepts SOE code 3 and DOE cartridges
- SOE code 7 and code B adaptors available
- Available in 316L stainless steel or Alloy 22 for maximum compatibility
- Wide range of connection size and types, cartridge lengths and materials of construction
- PED Compliant to Cat I
- Full traceability on all parts, material certification packs available
- Standard vessel rated to 28 bar g

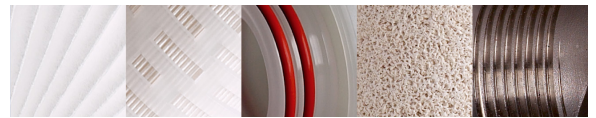

Accessories	Part Number
Pressure Gauge and Vent Valve Kit	53S010003
Ball Valve and Nipple	60S000010
Stand	50S010012
Mounting Bracket	50S010003
Spanner	53S000000

Features and Benefits

- Investment Cast head (avoids the porosity & inclusion problems associated with lower quality sand castings)
- Integral vent and drain sockets to aid efficient filter changeout complete with PTFE washer
- Unique interchangeable configuration enables choice of plug in or DOE cartridge styles in the same housing
- Every housing is etched with a serial number for full traceability
- Full housing customisation service available
- Available in both 316L Stainless Steel and Alloy 22

Product Range

Amazon Filters design and manufacture an extensive range of liquid filtration systems, providing one of the widest selections of filtration vessels, depth filters and pleated cartridges available today, ensuring excellent process and filter compatibility. As one of Europe's leading manufacturers of filtration systems, Amazon Filters is able to offer a solution to meet all filtration requirements. Amazons' filters are manufactured in clean room conditions and all products come with detailed Product Validation Guides and Technical Support documentation.



53A Series Technical Data

Maximum Operating Conditions

Pressure Range:	Standard
Liquid:	FV / 28 bar g
Gas Group 1:	FV / 20 bar g - 10" FV / 12 bar g - 20" FV / 9 bar g - 30"
Gas Group 2:	FV / 28 bar g
Temperature Range:	-10°C / 150°C

Materials of Construction

Body and Closure:	316L St.Steel Alloy 22
-------------------	---------------------------

Finish

Vessel:	Externally polished sump, bead blast head
---------	---

Connections

Inlet/Outlet:	1" BSPF
Vent:	¼" BSPF with Plug
Drain:	¼" BSPF with Plug

Model	Dimensions					Volume (Litres)	Dry Weight (kg)
	A	B	C	D	E		
10"	349	-	-	120	300	2.2	3.1
20"	603	-	-	120	554	3.7	4.1
30"	857	-	-	120	800	5.1	5.0

PED*/PESR** Status:

Compliant with PED 2014/68/EU and PESR 2016 No. 1105

Fluid Groups: Gas Group 1 + 2

Liquids Group 1 + 2

Marking: CE

Exceptions: Group 1 liquids - Explosive



ATEX:

Compliant with 2014/34/EU and

Explosive Atmosphere Regulations 2016

II 2G Ex h IIC T150 Gb

II 2D Ex h IIIC T150 Db

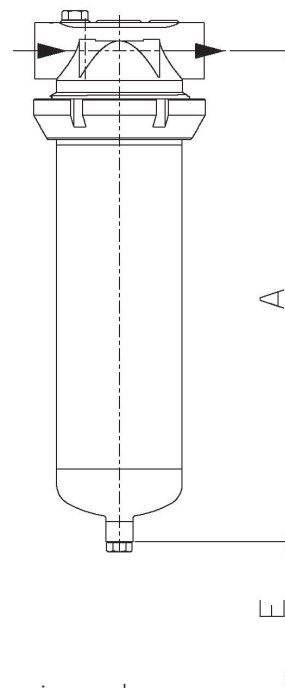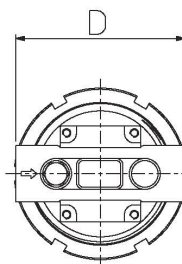

* Pressure Equipment Directive

** Pressure Equipment (Safety) Regulations No. 1105 Technical drawings for each individual housing showing dimensions, volume, weights, maximum temperatures and design pressures available on request.

Ordering Guide

53A1	1	0	X	N	10	B	1	C	V
No of Stacks	Cartridge Length	Cartridge End Cap Style	Closure Style	Mounting	Connection	Connection Type	Materials	Finish	Seal
53A1 - 1	0 - 5" 1 - Single 10" 2 - Double 20" 3 - Triple 30"	0 - Code 0 B - Code B 3 - Code 3 7 - Code 7 I - MCY1000 Insert*1	X - DIN Screwed Closure (Std)	N - None E - Bracket	10 - 1" 15 - 1 ½" 20 - 2"	B - BSPF (1" Only) C - EN1092-1 Flange PN16 E - EN1092-1 Flange PN40 H - ASME Flange 150# J - ASME Flange 300#	1 - 316 st.st (L Grade) 4 - Alloy 22	C - Externally Polished Sump, Bead Blast Head J - Externally Polished Sump, Bead Blast Head (Alloy 22 option only)	K - FFKM V - Viton S - Silicone

Example: 53A110XN10B1CV = 53 Series, 1 x 10" DOE cartridge, DIN Screwed closure, No mounting, 1" BSPF connections. St.Steel constructions with Viton gaskets, externally polished sump with bead blast head. *1 Can not be converted to accept code 3 cartridges.

AMAZON FILTERS LTD.

Albany Park Estate, Camberley, Surrey, GU16 7PG, ENGLAND

Tel: +44 (0) 1276 670 600 Email: sales@amazonfilters.co.uk Web: www.amazonfilters.com